

# Development and evaluation of a highly sensitive and specific radioimmunoassay for measurement of adrenocorticotropin levels in plasma

MANJULA K. GUPTA, PHD; TIM KOLAR, MS; LESLIE R. SHEELER, MD

■ The authors developed and evaluated a radioimmunoassay for adrenocorticotropin (ACTH) using Sep-Pak C<sub>18</sub> cartridges for the extraction of the hormone prior to assay. The assay has a sensitivity of 5 pg/mL. Inter- and intra-assay coefficients of variation were 9.1% and 4.8% at the level of 30 pg/mL (*n*=22). The mean ACTH level between 8 and 9 AM in 48 normal individuals was 20.3 ± 9.0 (S.D.) pg/mL with a range of 5–48 pg/mL. Mean (± S.D.) ACTH levels (8–9 AM) in nine normal subjects before and after 1 mg dexamethasone at midnight were 21 ± 8.8 and 9 ± 4.8 pg/mL, respectively, showing a mean suppression of 57% with a range of 38–76%. ACTH levels increased two- to fivefold in five normal subjects after a single oral dose of metyrapone. In 17 patients with hypercorticolism, ACTH levels indicated the status of the pituitary-adrenal axis.

□ INDEX TERMS: CORTICOTROPIN; RADIOIMMUNOASSAY □ CLEVE CLIN J MED 1988; 55:365–370

RECISE measurement of ACTH in plasma by radioimmunoassay (RIA) has been difficult, due primarily to the low plasma levels of this hormone and the nonspecific effects of the plasma matrix on most immunoassays. Although direct measurement of plasma ACTH by RIA has been described<sup>1-3</sup> and assays are currently available from a number of commercial enterprises,<sup>4-6</sup> most of these assays<sup>4-6</sup> lack sensitivity (lowest detectable dose > 25 pg/mL) and/or suffer from nonspecific interference caused by the plasma matrix.<sup>2</sup> Extraction of ACTH from plasma partially purifies the sample and allows it to be concentrated, and therefore increases assay accuracy and sensitivity. Various extraction procedures have been used,<sup>7-11</sup> of which

adsorption of the hormone to porous glass<sup>9</sup> or silicic acid<sup>10,11</sup> are the most common. Recently, the use of octadecylsilane-silica in Sep-Pak C<sub>18</sub> cartridges has been described for extraction of small polypeptide hormones, including vasopressin, angiotensin II, and ACTH.<sup>12-15</sup> We present results obtained with a plasma ACTH RIA that uses Sep-Pak C<sub>18</sub> cartridge extraction of the sample before measurement with a high-affinity ACTH antibody.

### MATERIALS AND METHODS

#### Materials

The antiserum used in this study was developed in a rabbit immunized with synthetic 1–24 ACTH coupled to bovine serum albumin. Purified porcine ACTH was used for iodination and standard. The chloramine T procedure of Hunter and Greenwood<sup>16</sup> was used for iodination. Briefly, 20 µL of 0.5 M phosphate buffer and 5 µL of

Departments of Immunopathology and Endocrinology, The Cleveland Clinic Foundation, Cleveland, Ohio. Submitted for publication June 1987; accepted Sep 1987.

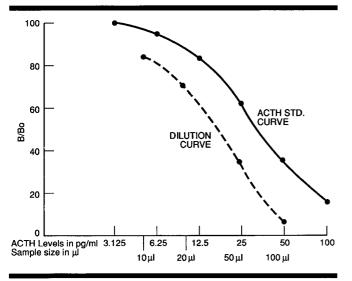


FIGURE 1. Inhibition of I-125-labeled ACTH bound to antibody by addition of unlabeled ACTH or plasma extract. Each point represents the mean of two determinations. The brackets indicate the S.D.

sodium I-125 (0.5 mCi [18.5 MBq]) were added to an iodination vial containing 1  $\mu$ g of ACTH (5  $\mu$ L) and gently mixed. Then 20  $\mu$ g of chloramine-T (10  $\mu$ L) was added and gently mixed for 15 seconds. The reaction was stopped by the addition of 20  $\mu$ L of sodium metabisulfite (15 µg). I-125-labeled ACTH was separated from unreacted I-125 on a 1 x 60-cm column of Sephadex G-50 Fine (Pharmacia, NJ), equilibrated, and developed in 0.06 M phosphate buffer, pH 7.4, containing 0.1% Triton X-100 (assay buffer) at 4°C. The flow rate was 10 mL/hr. Fractions of 1 mL were collected and monitored for radioactive counts. I-125-labeled ACTH was identified in fractions 28–31 and stored at –70°C until needed. For assay 100 µL of standard (3.0–100 pg/mL in assay buffer) or plasma extract was incubated with 100 µL of ACTH antibody diluted in assay buffer containing 1% normal rabbit serum (maximum binding 25-30%) for 24 hours at 4°C. One hundred microliters of freshly labeled or purified I-125-labeled ACTH was added to each tube and the incubation continued for 24 hours. Then 100 µL of rabbit anti-gamma globulin serum (diluted in assay buffer) was added and incubated for 4 hours. Antibody-bound and free ACTH were then separated by centrifugation at 5000 x g for 20 minutes at 4°C. The supernatants were immediately removed by aspiration and the precipitates counted in a gamma counter.

ACTH was extracted from plasma by use of Sep-Pak  $C_{18}$  cartridges (Waters, Milford, MA). The cartridge was

**TABLE 1**ASSAY CHARACTERISTICS

	_	
Sensitivity	=	5 pg/mL
$ED_{50}$ (dose required to reduce initial binding by 50%)	=	35 pg/mL
Precision:		
Interassay coefficient of variation		
at 7.1 pg/mL, $n = 22$	=	13.2%
at 28.8 pg/mL, $n = 22$	=	9.1%
Intraassay coefficient of variation		
at 7.1 pg/mL, $n = 4$	=	10.3%
at 28.8 pg/mL, $n = 4$	=	4.8%
Accuracy (recovery of known amount of ACTH		
added to plasma, $n = 22$ )	=	89 ± 5.3%

Defined as the lowest amount of unlabeled antigen tested that caused a significant inhibition of binding compared with zero dose; n = 10, p < 0.001.

first conditioned with 10 mL methanol then washed with 10 mL of distilled water. Then 0.5–1 mL of plasma was applied over 1 to 2 min and the cartridge washed with 10 mL of distilled water. ACTH was eluted with 6 mL of a 4:6 v/v mixture of TEAF buffer (1% formic acid titrated to pH 3.2 with triethylamine) and acetonitrile. The eluates were frozen in liquid nitrogen, lyophilized, and reconstituted in 0.5–1 mL of assay buffer just prior to assay. Plasma extracts, after lyophilization, can be stored indefinitely at  $-20^{\circ}$ C.

# Specimen collection

A 10-mL blood sample was drawn in a vacutainer containing 15 mg EDTA. Plasma was separated by centrifugation at 6900 x g at 4°C for 10 min and stored at -20°C until analyzed. Plasma from 48 normal subjects was obtained between 8 and 9 AM, and overnight dexamethasone suppression tests were performed on nine normal volunteers. For this test, a baseline sample was drawn between 8 and 9 AM for ACTH and cortisol assays, the subjects were given 1 mg of dexamethasone orally at midnight, and a sample for ACTH and cortisol levels was collected at the same time the following morning. In addition, single-dose metyrapone tests were performed in five normal subjects. Baseline samples for ACTH and 11deoxycortisol were drawn between 8 and 9 AM, the subjects were given 2 g of metyrapone orally at midnight, and another sample for ACTH and 11-deoxycortisol was drawn the following morning between 8 and 9 AM. Plasma cortisol and 11-deoxycortisol were measured by RIA.<sup>17,18</sup> ACTH was measured in samples collected from 17 patients with Cushing's syndrome (hypercorticolism); eight of these had pituitary Cushing's syndrome, four had adrenal tumor (three adenoma, one carcinoma), and five had ectopic ACTH-producing tumors.

#### RADIOIMMUNOASSAY OF ADRENOCORTICOTROPIN ■ GUPTA AND ASSOCIATES

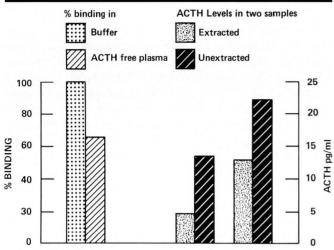


FIGURE 2. Effect of plasma matrix on antibody binding. I-125-labeled ACTH bound to antibody in buffer (0.06 M phosphate buffer, pH 7.4) medium is taken as 100% (left scale). On the right scale, ACTH levels as detected in two samples before and after extraction.

#### RESULTS

Figure 1 shows the inhibition of I-125-labeled-ACTH binding by added unlabeled ACTH, and Table 1 lists the analytical parameters for this assay. The sensitivity of the assay is 5 pg/mL when a 0.5-mL aliquot of plasma is used for extraction. The sensitivity can be further enhanced by use of 1 to 2 mL plasma, thereby concentrating the hormone two- to fourfold without loss of extraction efficiency. Dilutions of plasma extract showed parallel inhibition of binding with the ACTH standard. The mean ACTH level in 48 normal subjects was  $20.3 \pm 8.9$ pg/mL (S.D.) with a range of 5 to 48 pg/mL. The mean recovery of ACTH added to 1 mL plasma (n = 22) at the level of 25 pg/mL was  $89 \pm 5.4\%$  (S.D.). The interassay coefficient of variation (n = 22) at the level of 7.1 pg/mL was 13.2% and at 28.8 pg/mL was 9.1%. The corresponding intra-assay coefficients of variation were 10.3% and 4.8%.

*Figure 2* illustrates the effect of plasma matrix on antibody binding. One hundred microliters of ACTH-free plasma from a normal subject (postdexamethasone, ACTH removed by Sep-Pak) caused 35% inhibition in antibody binding. This interference was confirmed by running ten samples before and after extraction in the assay. Significantly higher values (p < 0.0001) were obtained in the unextracted samples. The mean of unex-

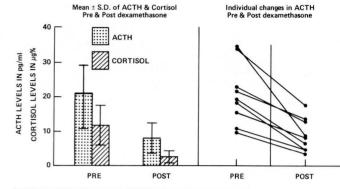


FIGURE 3. ACTH and cortisol levels in normal subjects pre- and post-dexamethasone. Individual ACTH values in these subjects are shown on the right.

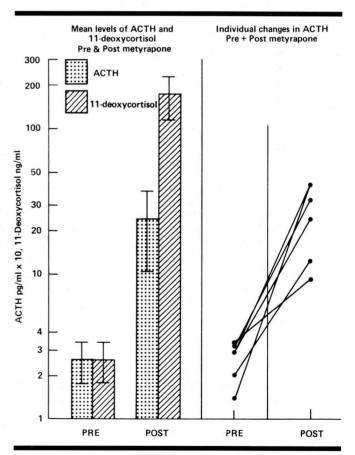


FIGURE 4. ACTH and 11-deoxycortisol post-metyrapone in five normal subjects. Mean  $\pm$  S.D. are shown on the left side and individual changes in ACTH are shown on the right side.

## RADIOIMMUNOASSAY OF ADRENOCORTICOTROPIN ■ GUPTA AND ASSOCIATES

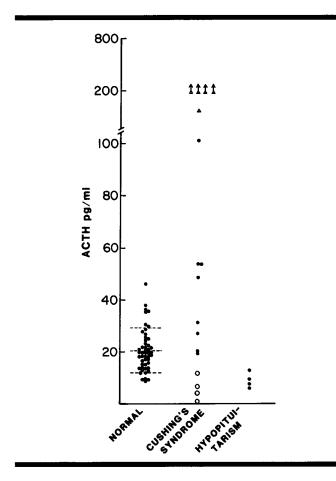


FIGURE 5. ACTH levels in 48 normal subjects (the broken lines represent mean and 2 S.D. range), 17 patients with Cushing's syndrome (adrenal tumor [open circles]), excess production of ACTH by the pituitary (closed circles), and ectopic ACTH production (triangles), and four patients with hypopituitarism.

tracted samples was  $31 \pm 8.8\%$  (S.D.) and mean for extracted samples was  $14.5 \pm 5.4\%$  (S.D.). The degree of increase varied between 19% and 380% from sample to sample, indicating the need of extraction for meaningful information.

Figure 3 shows the plasma ACTH and cortisol levels in normal subjects before and after dexamethasone. A postdexamethasone decrease in ACTH, which paralleled cortisol levels, was observed in all nine subjects. The mean basal ACTH level was 20.7  $\pm$  8.3 pg/mL (S.D.), postdexamethasone was 8.8  $\pm$  4.6 pg/mL (S.D.), and percent suppression was 57  $\pm$  13% (S.D.). The corresponding cortisol levels were 12.4  $\pm$  5.8 µg/dL (S.D.) and < 2.5 µg/dL, respectively, representing > 80% suppression.

Figure 4 shows the ACTH and 11-deoxycortisol levels in five normal subjects during a single-dose metyrapone test. The mean ( $\pm$  S.D.) baseline values were 25.7  $\pm$  8.3 pg/mL for ACTH and  $2.5 \pm 0.8$  ng/mL for 11-deoxycortisol, which increased to  $237 \pm 139$  (S.D.) pg/mL and 172 ± 56 ng/mL (S.D.), respectively, after metyrapone. Figure 5 shows ACTH levels in 48 normal subjects and 17 patients with hypercorticolism as documented clinically by cortisol determinations. Of these 17 patients, eight patients had pituitary Cushing's disease, four patients had adrenal adenoma or cancer, and five patients had ectopic ACTH-producing tumors. The ACTH levels in eight patients with pituitary Cushing's syndrome ranged between 20 and 108 pg/mL with the mean and S.D. of  $45.3 \pm 29$  pg/mL. In four patients with adrenal adenoma (one with cancer), the mean  $(\pm S.D.)$  ACTH level was  $7.2 \pm 3.4$  pg/mL with a range of < 5 to 12 pg/mL, and in five patients with ectopic ACTH-producing tumors, the mean ( $\pm$  S.D.) ACTH level was 500  $\pm$  406 pg/mL with a range of 160 pg/mL to 11,000 pg/mL. In three patients with hypopituitarism, the ACTH level averaged  $10.6 \pm$ 2.8 pg/mL (S.D.). In one patient with Nelson's syndrome, the ACTH level was 900 pg/mL.

### DISCUSSION

Measuring ACTH in biological fluids is important in assessing hypothalmopituitary adrenal function and in the differential diagnosis of hypercorticolism. To do this both normal and subnormal ACTH levels must be detected, which requires a highly sensitive and specific assay. The initial commercial ACTH assays suffered from lack of sensitivity and required preliminary extraction from 5 to 10 mL plasma to concentrate ACTH in order to detect normal levels.19 They were also tedious and lacked precision. In recent years, a number of direct assays have been made available by commercial enterprises.<sup>4-6</sup> These assays, although simple to perform, lack sensitivity to detect subnormal and normal levels and have a variable degree of specificity. False-high values (ACTH as high as 300 or 700 pg/mL) were detected in normal individuals by use of one of these assays<sup>4</sup> in our laboratory. Similar false-high levels have been reported by others.<sup>20</sup> The false-high levels could be due to plasma peptidases, which may degrade radioiodinated ACTH, or to other substances in plasma capable of inhibiting antigen-antibody binding. More recently, Nicholson et al<sup>2</sup> reported a rapid RIA procedure in unextracted human plasma. They also observed nonspecific inhibition of antigen-antibody binding by plasma components, contributing to the apparent immunoreactive ACTH content of the samples. But, unlike us, they observed a relatively constant amount of inhibition by most of the plasma samples. To compensate for this nonspecific effect, they used carefully selected lots of bovine serum albumin, which were extracted with silicic acid to remove ACTH-like activity.

We found that the ability of ACTH-free plasma to inhibit antibody binding varied significantly from sample to sample and could not be compensated by any one lot of bovine serum albumin. These data suggested the need for extraction of ACTH from plasma. Various methods of ACTH extraction from plasma have been applied<sup>7-11</sup> and those most commonly used are based on adsorption of the ACTH to porous glass or silicic acid.9-11 Recently, octadecylsilane-silica Sep-Pak in C<sub>18</sub> cartridges has been used for the separation of vasopressin, angiotensin II, and ACTH from plasma proteins prior to RIA.<sup>12-15</sup> We used these cartridges for the separation of ACTH from plasma and found the normal range to be in accordance with previously reported values employing other extraction methods<sup>19</sup> but significantly lower than the range obtained by direct measurement of plasma ACTH in our lab or results reported with other commercial reagents.<sup>20</sup> These higher normal values in the direct assay also suggest nonspecific interference.

So far our data on normal individuals and patients with various diseases show that the ACTH levels as

#### REFERENCES

- Berson SA, Yalow RS. Radioimmunoassay of ACTH in plasma. J Clin Invest 1968; 47:2725–2751.
- Nicholson WE, Davis DR, Sherrell BJ, Orth DN. Rapid radioimmunoassay for corticotropin in unextracted human plasma. Clin Chem 1984; 30:259–265.
- Kao PC, Jiang N-S, Carpenter PC. Human corticotropin (ACTH) radioimmunoassay with synthetic 1–24 ACTH. Clin Chem 1979; 25:1267–1273.
- ACTH <sup>125</sup>I radioimmunoassay kit: Direct radioimmunoassay for ACTH in plasma. Package Insert, Immuno Nuclear Corp, Jan 1984.
- ACTH<sup>125</sup>I radioimmunoassay kit: Direct radioimmunoassay for ACTH in plasma. Package Insert, Diagnostic Systems Laboratories, Jan 1984.
- ACTH <sup>125</sup>I radioimmunoassay kit: Direct radioimmunoassay for ACTH in plasma. Package Insert, CIS-Sorin (distributed by Damon Diagnostics, Maine 02194), 1983.
- Demura H, West CD, Nugent CA, et al. A sensitive radioimmunoassay for plasma ACTH levels. J Clin Endocrinol Metabol 1966; 26:1297–1302.
- Orth DN, Island DP, Nicholson WE, et al. ACTH radioimmunoassay: Interpretation, comparison with bioassay, and clinical application. [In] Hayes RC, Goswitz FA, Murphy BEP, eds. Radioisotopes in Medicine: In Vitro Studies. CONF-671111. Springfield IL, Clearinghouse for Federal Scientific and Technical Information, National Bureau of Standards, US Dept of Commerce, 1968, pp 251–272.
- 9. Ratcliffe JG, Edwards CRW. The extraction of adrenocorticotropin

measured by this assay generally reflect accurately the status of the subject's pituitary-adrenal axis. Our results obtained during treatment with dexamethasone show that this assay can detect significant suppression of ACTH in normal subjects. This may provide more accurate assessment of ACTH status in depression than cortisol levels. Cortisol levels have been used conventionally to diagnose hypercorticolism but the differential diagnosis of this condition is very difficult even with extensive testing. ACTH levels measured in 17 patients with Cushing's syndrome could be used to separate cases where hypercortisolism was caused by excess ACTH production by the pituitary from those whose Cushingnoid state was due to ectopic ACTH production; however, they could not be used to identify all patients with Cushing's syndrome secondary to an adrenal tumor as having suppressed plasma ACTH levels. Thus, although the number of patients tested was small, the data suggest that plasma ACTH assays as performed in our laboratory may aid in differential diagnosis of hypercorticolism.

> MANJULA K. GUPTA, PHD Department of Immunopathology The Cleveland Clinic Foundation One Clinic Center 9500 Euclid Avenue Cleveland, Ohio 44195

and arginine-vasopressin from human plasma by porous glass. [In] Kirkham KE, Hunter WM, eds, Radioimmunoassay Methods. Ed-inburgh, Scotland, Churchill Livingstone, 1971, pp 502–512.

- Kreiger DT, Allen W. Relationship of bioassayable and immunoassayable plasma ACTH and cortisol concentrations in normal subjects and in patients with Cushing's disease. J Clin Endocrinol Metabol 1975; 40:675–687.
- Orth DN. Adrenocorticotropic hormone (ACTH). [In] Jaffe BM, Behrman HR, eds. Methods of Hormone Radioimmunoassay. 2nd ed. New York, Academic Press, 1979, pp 245–283.
- 12. Chrousos GP, Heinrich M, Schulte M, et al. The corticotropinreleasing factor stimulation test. An aid in the evaluation of patients with Cushing's syndrome. New Engl J Med 1984; 310:622–626.
- Ysewijn-Van Brussel KARN, De Leenheer AP. Development and evaluation of a radioimmunoassay for Arg8-vasopressin, after extraction with Sep-Pak C<sub>18</sub>. Clin Chem 1985; 31:861–863.
- Suzuki H, Ferrario CM, Speth RC, Brosnihan KB, Smeby RR, DeSilva P. Alterations in plasma and cerebrospinal fluid norepinephrine and angiotensin II during the development of renal hypertension in conscious dogs. Hypertension 1983; 5:I-139–I-148.
- Oldfield EH, Schulte HM, Chrousos GP, et al. Corticotropin-releasing hormone (CRH) stimulation in Nelson's syndrome: Response of adrenocorticotropin secretion to pulse injection and continuous infusion of CRH. J Clin Endocrinol Metabol 1986; 62:1020–1026.
- Hunter WM, Greenwood FC. Preparation of iodine-131 labeled human growth hormone of high specific activity. Nature 1962; 194:495–496.

## RADIOIMMUN0ASSAY OF ADRENOCORTICOTROPIN ■ GUPTA AND ASSOCIATES

- Brock P, Eldred EW, Woiszwillo JE, Doran M, Schoemaker HJ. Direct solid-phase <sup>125</sup>I radioimmunoassay of serum cortisol. Clin Chem 1978; 24:1595–1598.
- Perry LA, Al-Dujaili EAS, Edwards CRW. A direct radioimmunoassay for 11-deoxycortisol. S Thyroids 1982; 39:115–128.
- ACTH radioimmunoassay kit (code IM.66) Package Insert. Amersham/ Searle Corp., Arlington Hts., IL 60005.
- Hedner P, Norden NE, Valdemarsson S. Falsely high values for corticotropin by the CIS-Sorin RIA method. Clin Chem 1981; 27:1146–1147.

